

**Nicolet**  
**StethoDop™**  
VASCULAR DOPPLER



- \* **Twists on and off most stethoscopes**
- \* **Convenient, lightweight & versatile**
- \* **Crystal-clear signal**

**VIASYS™**  
HEALTHCARE  
*Excellence For Life*  
NEUROCARE GROUP

# Nicolet StethoDop™

## Turn your Stethoscope into a Doppler

Your most familiar diagnostic tool just got better. Now Nicolet Vascular offers something truly unique for your stethoscope - a 1 .5 ounce quick-connect Doppler. Nothing is more convenient, lightweight, powerful and versatile than the StethoDop. It is the first Doppler that twists on and off most bell-style stethoscopes in seconds. Most importantly, it never interferes with the normal operation of your stethoscope.

### Compact Convenience

The ingenious StethoDop provides many benefits previously available only in larger Doppler systems, including crystal-clear signals, high-quality audio and a sensitive, easy-to-use 5 MHz transducer. The pre-angled, wide beam crystals cover a large area, allowing for faster vessel location. Simply align the StethoDop probe face flat along the desired vessel, press the instant “on” button and listen to something new.

### Powerful diagnostic aid


Use the StethoDop to quickly and easily listen to blood flow throughout the peripheral vascular system.

- Obtain systolic blood pressures even with weak pulses.
- Perform fast ABI (ankle/brachial index) exams to screen for peripheral arterial disease.
- Evaluate blood flow to wound sites.
- Quickly screen for lower extremity vascular conditions.
- Use as an aid to catheter insertion and to determine vessel patency.

Now sensitive blood flow detection is available in even the noisiest situations, including emergency vehicles, aircraft, trauma centers and other high-traffic areas.

### Exceptional Ownership Experience

VIASYS NeuroCare captures the knowledge and expertise of four companies (Nicolet Vascular, Nicolet Biomedical, Grason-Stadler and Toennies) representing over five decades of ingenuity and innovation. Nicolet Vascular is a leading manufacturer of hand-held and desktop Dopplers for obstetric and vascular applications. We strive to build upon our hard-earned reputation and continue to meet the rigorous standards for quality and reliability that our customers have come to expect. Nicolet StethoDop follows in the VIASYS tradition of technological innovation and superior customer service. The StethoDop is manufactured in our ISO 9001 certified facility and is covered by a 5-year parts and 1-year labor warranty.

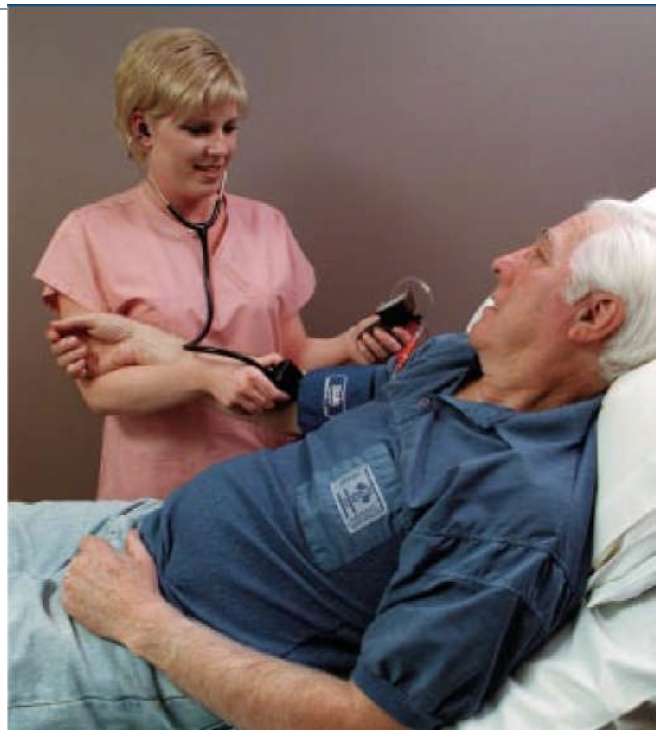

### Catalog Numbers

StethoDop .....SD50  
StethoDop with cardiology-style stethoscope.....SD52  
**StethoDop Accessories:**

Battery/Non-chill ring/speaker coupler pack .....A600  
Carrying Pouch .....A410  
Cardiology-style stethoscope ..... SS1 0

### Technical Specifications

Weight: ..... 45 grams (1.5 oz.)  
Dimensions: .....6.25 x 3.75 x 2.5 cm. (2.5 x 1.4 x 1 in.)  
Battery Type: .....3 each, 3 VDC lithium-type 1/3N cells



Distributed By:

GLOBAL MEDICAL DISTRIBUTION LTD

[www.gmdltd.com](http://www.gmdltd.com)